


Obstruction Removal Kit for ST Nebulizer Nebulizers

On occasion the ST nebulizers become clogged. This kit is specifically designed for PFA-ST and PolyPro-ST nebulizers.

- ▲ Includes a guide tool and 10 pieces of fused silica.
- ▲ 0-3 ml/min sample flow rate with peristaltic pumping.
- ▲ Easy to use:
 - tighten guide onto back of ST nebulizer
 - slide fused silica through.
- ▲ Silica refills also available.

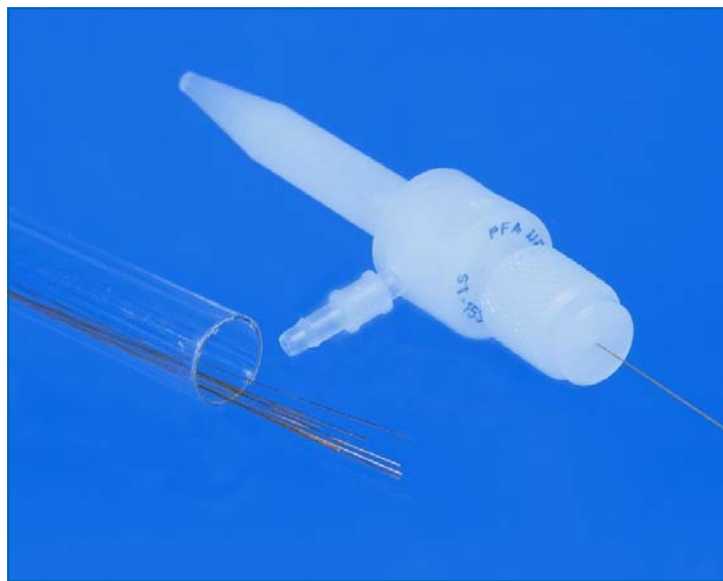

Obstruction Removal Kit shown with

