


## PFA High Purity Sample Introduction System for GV ICP-MS Systems

Inert sample introduction technology for improved background and detection limits

An easy to install, inert, high purity, o-ring free upgrade kit which includes two PFA MicroFlow nebulizers, PureCap endcap, PureChamber PFA spray chamber, adapter, o-ring free injector and precision quartz torch. A reliable ultra pure analysis technology to complement GV ICP-MS systems

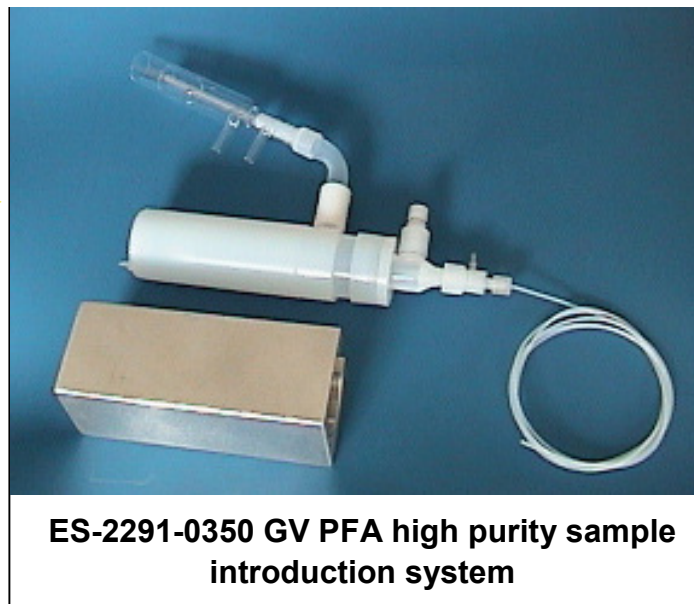

**ES-2291-0350 GV PFA high purity sample introduction system**

- Easy installation on all GV ICP-MS systems.
- Clean Teflon® PFA components for lower background equivalent concentration (BEC) and detection limits.
- O-ring free connections —ensures lowest BEC's, highest chemical resistance and easy cleaning.
- Chemically resistance—suitable for nearly all samples, including strong acids, alkalis and organic solvents.
- Sapphire or platinum injector.
- Two self-aspirated or pumped PFA nebulizers included.
- Low sample consumption—ideal for VPD and preconcentrated samples
- Optional Aluminum cooling block for improved stability and low maintenance.

Elemental Scientific Inc.  
2440 Cuming St.  
Omaha, NE 68131 USA  
[www.elementalscientific.com](http://www.elementalscientific.com)

Phone: 402.991.7800  
Fax: 402.991.7799  
Email: [esi@icpms.com](mailto:esi@icpms.com)



